


How to create the Patients.txt file for Lab Importer Pro manually in Praxis


Download the Praxis Query file and unzip contents.

<http://www.schmitzsoftware.com/LabImporterPro/Instructions/PraxisQuery.zip>


Step 1. Log into Praxis and open the Query tab – If you have the Create_Patients.txt_File.p3q query installed, go to Step 4

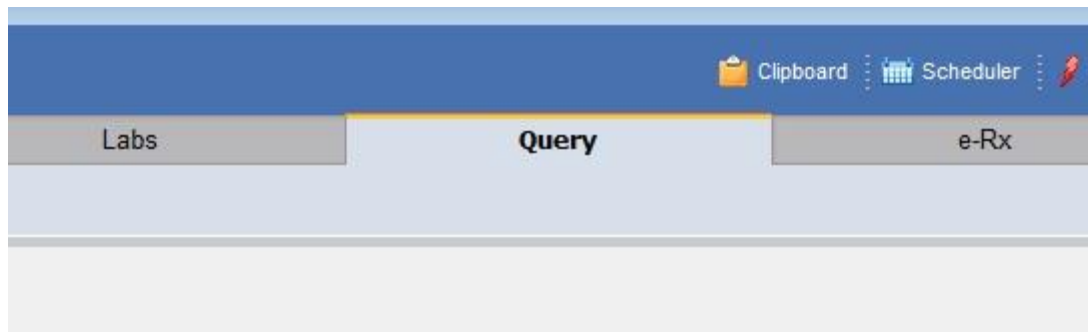


Step 2. Click Import

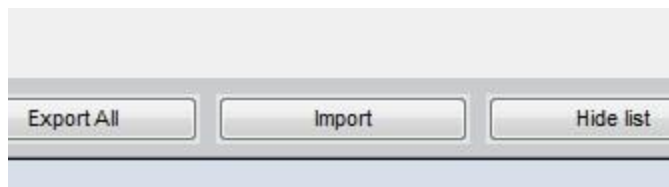


Step 3. Select the Create_Patients.txt_File.p3q from the download location.

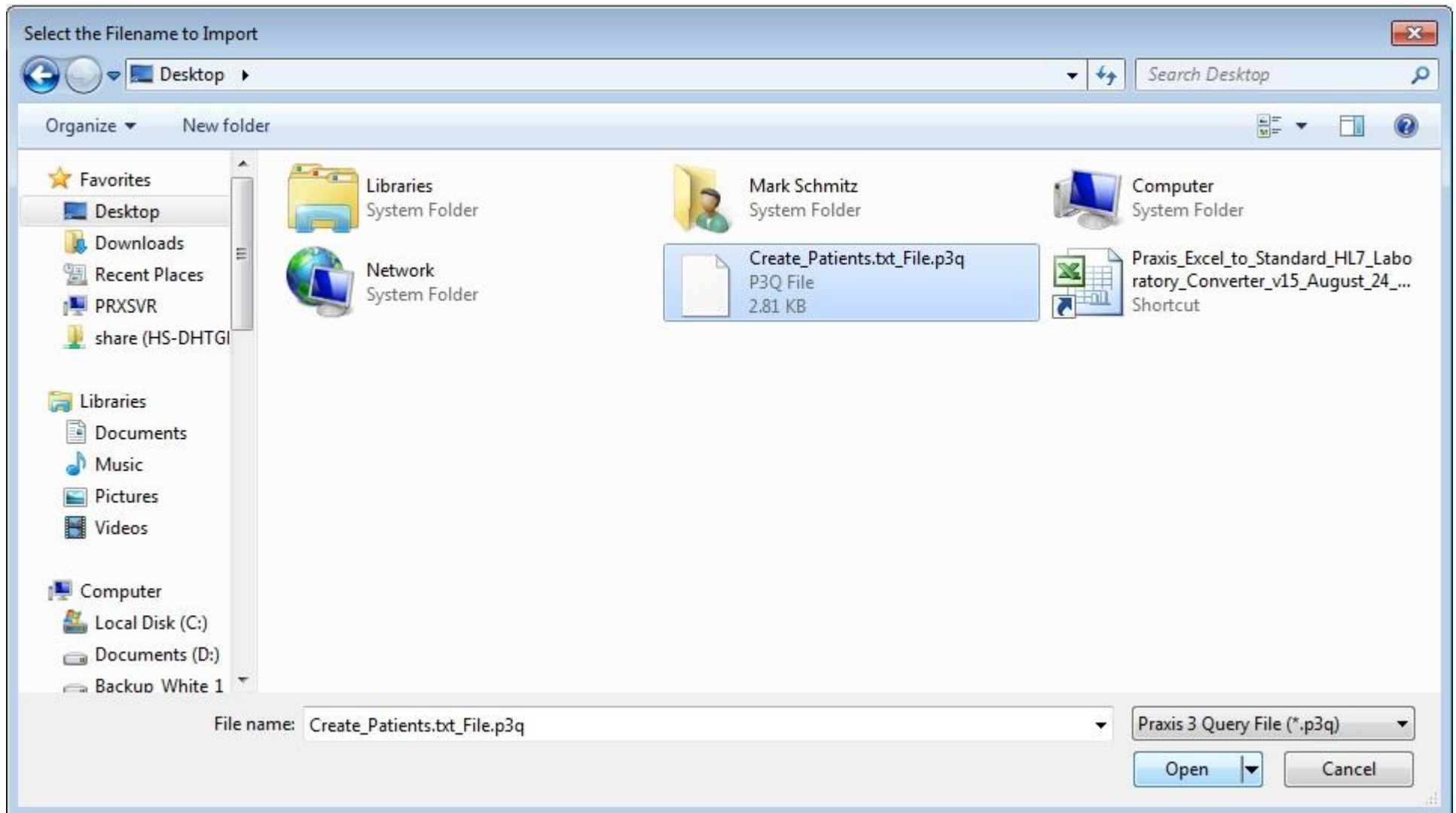


You should get the following message:

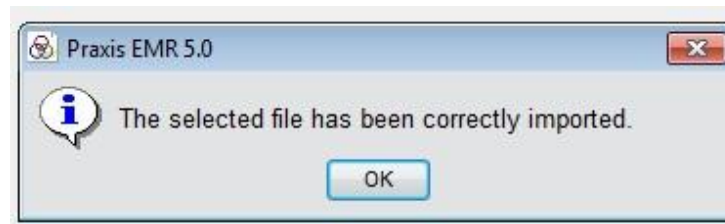


Step 4. Search for the Create Patients.txt File query

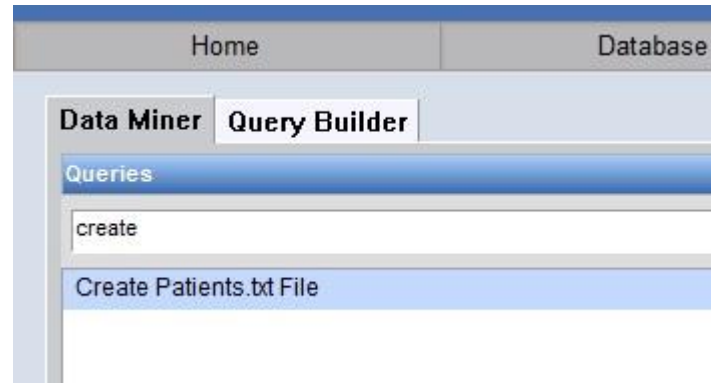


Step 5. Double click on the Create Patients.txt File query to load

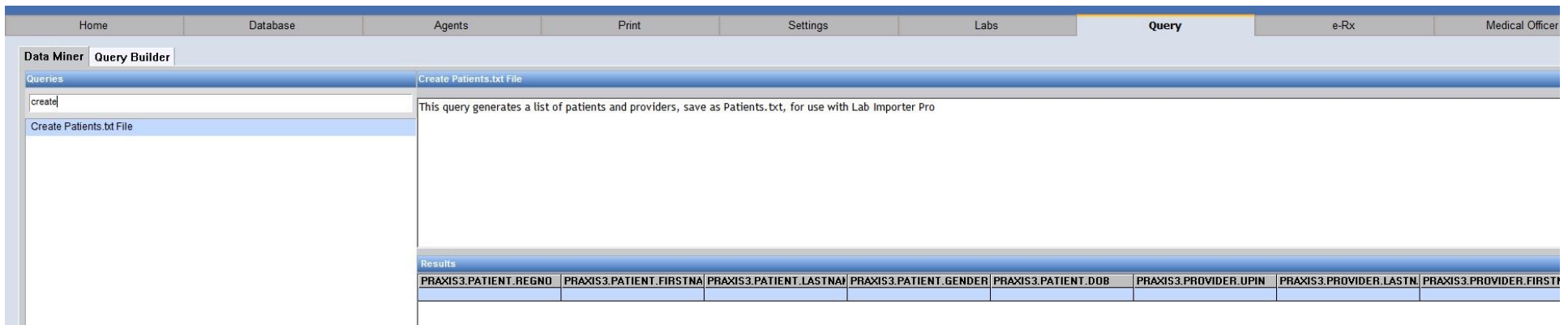


Step 6. Click Execute

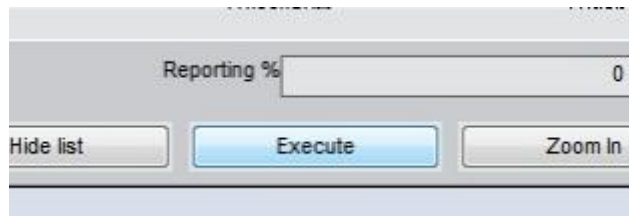


Step 7. Right click on the results and select Export data

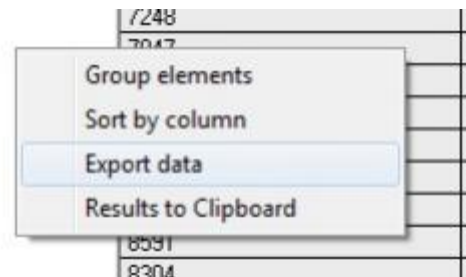

Step 8. Go to the directory that you have the program LabImporterPro.exe and save the file as Patients.txt

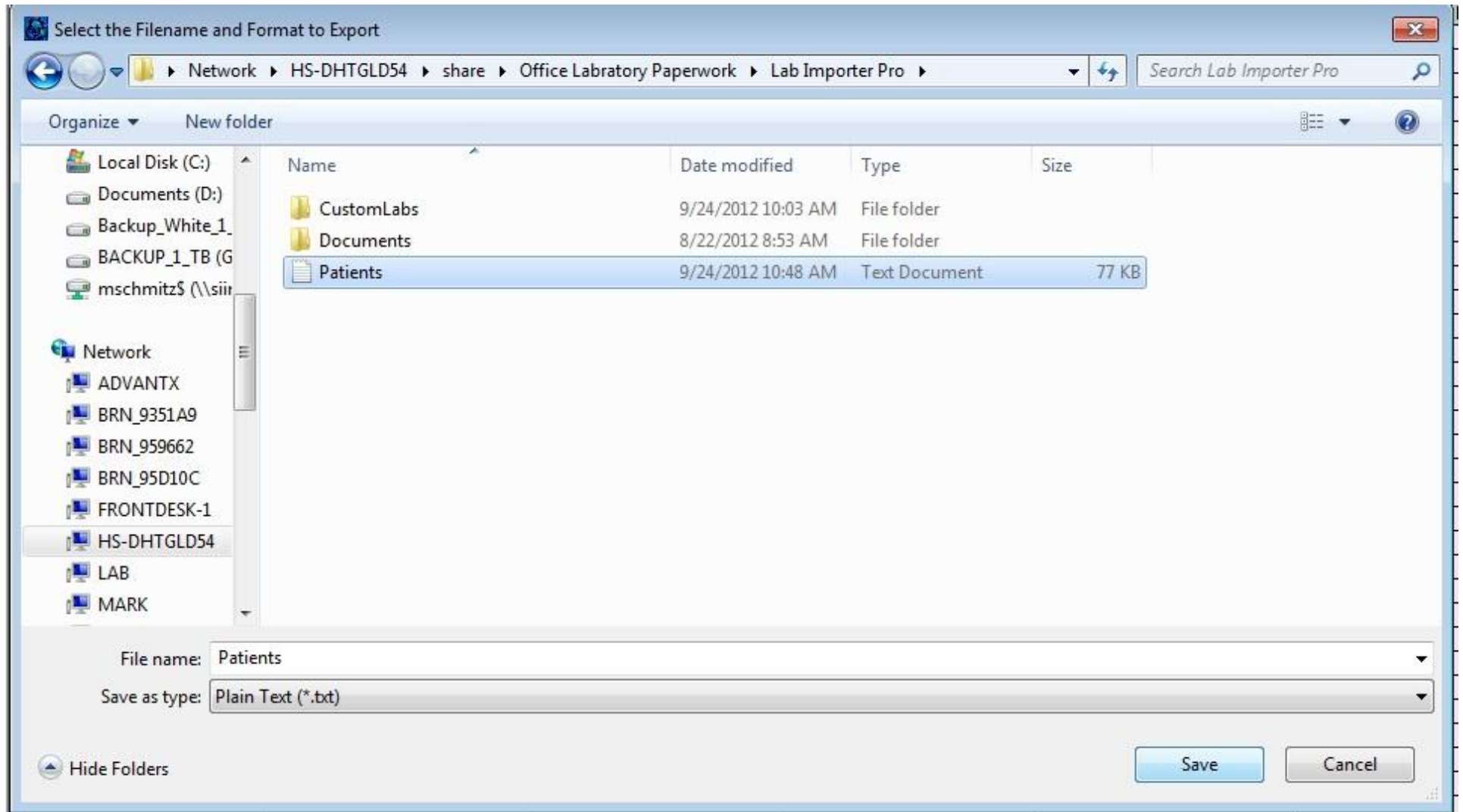

Step 9. If you get the following message, click Yes

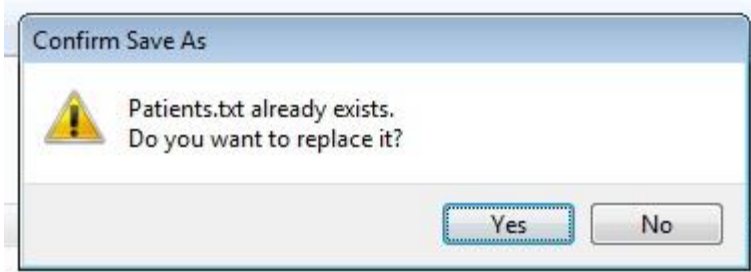

Once you see the following message, click OK. Your Patients.txt file is updated.

